

COTTON PRODUCTS



To meet your need for perfection

... all things cotton

roeko

- Reliable high manufacturing standard – qualityfocussed and fully automated
- Pharmaceutical quality raw materials – tested for purity, foreign fibres, absorbency and surfactant residues
- Versatile and well thought out – adapted to the needs of dentists

Cotton products Made in Germany

roeko – one name – long-standing experience – extensive technical know-how – utmost attention to quality

The quality of cotton products definitely pays dividends if one does the sums correctly.

The reliable high manufacturing standard of ROEKO cotton products contributes to lower consumption,

trouble-free working and hence to truly perfect results.



LUNA GOLD

See and feel the difference

LUNA GOLD is encased in skin-friendly mesh and can be removed without prior wetting – this prevents mucosal reaction and saves time. The gold stripe clearly differentiates LUNA GOLD from conventional dental rolls.

For regular use we recommend

LUNA – long-lasting absorbency and exceptional dimensional stability

COTTON-ROLL – particularly soft with rapid absorption.

ROEKO dental rolls are bleached chlorine-free.

ROLL-O-MAT

The system is both simple and ingenious: slight pressure on the lid is sufficient to dispense a single cotton roll and keep it positioned reliably on the dispenser tip. The ROLL-O-MAT is available as plastic dispenser or in a stainless steel version.



LUNA FOR KIDS

Specially designed for the very small patients

LUNA FOR KIDS dental rolls are particularly short and yet dimensionally stable like their larger Luna counterparts. This leaves plenty of room to work, even in the smallest of mouths. ROEKO dental rolls are bleached chlorine-free.


LUNAMAT

The LUNAMAT is equipped with a drawer which dispenses 2-3 cotton rolls when pulled out.



DENTAL ROLL DISPENSERS

ROEKO dental roll dispensers keep cotton rolls within reach during treatment while at the same time protecting them hygienically from spray and dust.

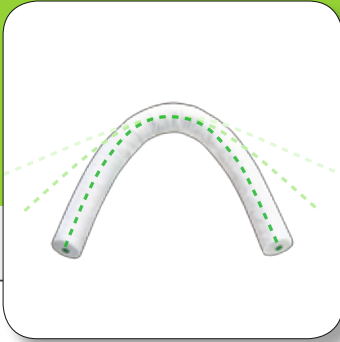

PAROTISROLL

Multipurpose cotton roll with spring action

The spring rod inside the cotton roll provides a spring action and makes the PAROTISROLL self-fixating. When used buccally, lingually or labially, it keeps the mouth open, holds down the tongue and retracts the lips. When positioned appropriately, it also prevents the flow of saliva through pressure on the parotid gland.

PAROTISROLL ANATOMIC

By Dr. Hannelore Oberschachtsiek PAROTISROLL anatomic is the Parotisroll with flattened ends. The flattened ends facilitate placement of the cotton roll and reduce the foreign body feeling in the mouth.

PAROTIMAT

For convenient individual dispensing and hygienic storage of 50 PAROTISROLLS.



COTTON PELLETS

Uniform in shape and size


The fully automated production of the ROEKO cotton pellets guarantees their uniform shape and size and thus, reliability. They are the ideal aid for drying cavities or applying medications.

#000	#00	#0	#1	#2
○	○	○	○	○
Ø 3 mm	Ø 4 mm	Ø 4,8 mm	Ø 6,3 mm	Ø 7,6 mm

SOLOMAT

The elegant dispenser

In this stainless steel dispenser, the cotton pellets are spring-loaded to the dispensing surface. The scratch- and abrasion-resistant surface prevents bacteria from adhering and is easy to clean and disinfect.

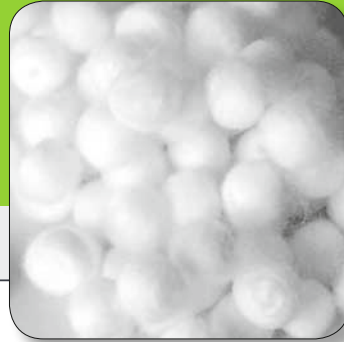

SOLOMAT-N

For individual dispensing of cotton pellets or foam pellets

Autoclavable plastic dispenser with refill system. The cotton pellets are moved manually to the dispensing surface.



COMPRECAP

COMPRECAP compression caps are a simple and intelligent aid for preparing impression taking. They stop bleeding – by compression and dry tooth stumps.

COMPRECAP ANATOMIC

COMPRECAP ANATOMIC are ideally adapted to the anatomy of the dental arch through semi-circular cut-outs.

SOLOMAT-CLICK

The smart dispensing system

SOLOMAT-CLICK consists of a dispensing grille, cover and refill. These 3 parts are clicked together when repurchasing a refill. The stainless steel base gives the dispenser the necessary stability.



COTTON PELLETS DISPENSER


The ROEKO cotton pellet dispensers keep cotton pellets within reach during treatment while at the same time protecting them hygienically from spray and dust. A firm base and well-designed lid allow dispensing cotton pellets with only one hand. The rounded rods of the dispensing area prevent the pellets from being squeezed on removal as opposed to a dispensing area with a grid structure.

Order information

	Contents	REF
LUNA	300 g	
	# 1	100 001
	# 2	100 002
	# 3	100 003
	# 4	100 004
	500 pcs.*	
	# 1	105 001
	# 2	105 002
	# 3	105 003
	1000 pcs.*	
	# 1	104 001
	# 2	104 002
	# 3	104 003
LUNA FOR KIDS	200 pcs. # 2	100 032
LUNA GOLD	300 g	
	# 1	109 121
	# 2	109 122
	# 3	109 123
	1000 pcs.*	
	# 1 # 2 # 3	109 101 109 102 109 103
COTTON-ROLL	1000 pcs.* # 1 # 2 # 3	104 011 104 012 104 013
ROLL-O-MAT plastic, white 134 × 100 × 48 mm	1 pc.	120 004
	chrome nickel steel 134 × 100 × 43,5 mm (can be used with cotton roll sizes 1–3), standard length (36 – 40 mm). Not suitable for LUNA FOR KIDS.	120 005
LUNAMAT plastic, white 52 × 106 × 100 mm (can be used with cotton roll sizes 1–4), standard length (36 – 40 mm) and with LUNA FOR KIDS.	1 pc.	120 000

* Only for export

© 2016 Coltene/Whaledent AG - www.coltene.com



			Contents	REF
PAROTISROLL	ø	L		
	# 1	9 mm	80 mm	100 pcs. 160 001
	# 2	10 mm	80 mm	160 002
	# 3	10 mm	100 mm	160 003
	# 5	8 mm	100 mm	160 005
	anatomic	10 mm	100 mm	50 pcs. 160 013
PAROTIMAT 112 × 103 × 82 mm			1 pc.	162 100
COTTON PELLETS	ø			
	3.0 mm	# 000	1.5 g	210 000
	4.0 mm	# 00	4.0 g	210 100
	4.8 mm	# 0	4.0 g	210 110
	6.3 mm	# 1	10.0 g	210 111
	7.6 mm	# 2	10.0 g	210 112
SOLOMAT¹ stainless steel 72 × 53 × 50 mm		# 000	1 pc.	6001 9711
		# 00		6001 9712
		# 0		6001 9713
		# 1		6001 9714
		# 2		6001 9715
SOLOMAT-CLICK² plastic 43 × 50 × 41 mm		# 000	1 pc.	224 000
		# 00		224 100
		# 0		224 110
		# 1		224 111
REFILLS		# 000	3 pcs.	214 000
		# 00		214 100
		# 0		214 110
		# 1		214 111
SOLOMAT-N² plastic 75 × 43 × 52 mm		# 000	1 pc.	223 000
		# 00		223 100
		# 0		223 110
		# 1		223 111
		# 2		223 112
REFILLS		# 000	6 pcs.	213 000
		# 00		213 100
		# 0		213 110
		# 1		213 111
		# 2		213 112
COMPRECAP	# 1	120 pcs.		530 008
	# 1	600 pcs.		531 008
	# 3	120 pcs.		530 010
	# 3	600 pcs.		531 010
	# 5	60 pcs.		530 014
	# 5	400 pcs.		531 014
	Set (#1/3/5)	120 pcs.		530 000
COMPRECAP ANATOMIC	# 1	120 pcs.		531 201
	# 3	120 pcs.		531 203
	# 5	60 pcs.		531 205
	Set (#1/3/5)	120 pcs.		531 299

¹⁾ Dispenser is supplied without cotton pellets

²⁾ Dispenser is supplied with cotton pellets

Coltene/Whaledent Ltd.
The President Suite
Kendal House, Victoria Way
Burgess Hill, West Sussex
RH15 9NF / UK
T +44 1444 235486
F +44 1444 870640
info.uk@coltene.com

Coltene/Whaledent AG
Feldwiesenstrasse 20
9450 Altstätten/Suisse
T +41 (0)71 757 53 00
F +41 (0)71 757 53 01
office@coltenewhaledent.ch

Coltene/Whaledent Inc.
235 Ascot Parkway
Cuyahoga Falls, Ohio 44223/USA
T +1 330 916 8800
F +1 330 916 7077
productinfo@coltenewhaledent.com

